

User Guide

Reference Ruler


User Level Radiologists

Overview

There is now a tool where users can add a reference ruler into the diagnostic viewer, the same thing seen on other PACS. Regardless of the image, the viewer will have a reference ruler on the side and bottom.

How to Use

On the diagnostic viewer, look for the Reference ruler tool on the tools list. Once added, click to enable/disable the reference ruler

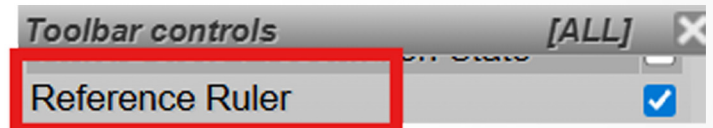

NOTE: Accurate reference measurements can only be generated when the DICOM Pixel Spacing tag is available.

Additional Resources

Click the button below to watch a short video tutorial demonstrating how to use this feature:

[Click Here](#)